

710

Exam Chair Doubles as a
Minor Procedure Table

800.735.0357
haag-streit-usa.com

Pillow Headrest

Backlit Membrane Switches

Reclines 180° for Minor Procedures

Ergonomic Design Speeds Up Exam


Power Base and Articulating Top



710

Exam/Procedure Chair

800.735.0357
haag-streit-usa.com

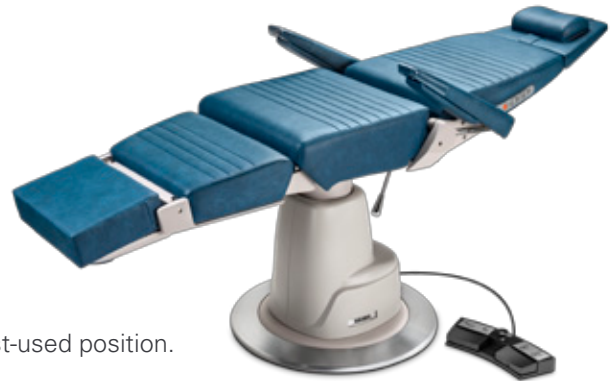

● Ergonomic Features

Round base optimizes access to patient.

Memory programming lets you pre-select the chair's most-used position.

Dual-rotation locks allow chair to swivel 330° from both sides of the chair.

Large hinged footrest safely supports patient's full weight.

Padded armrests lock in place and release with a single hand.

Comfort construction reduces patient stress and shortens procedure times.

Detached footswitch allows for hands-free control of all chair movements.

● Pillow Headrest

Standard non-articulating pillow headrest. Optional articulating version features a sliding, removable magnetic pad and fold-away design for otology – exclusive to the 710 chair.

● Power Base and Articulating Top

Finger switches and a corded footswitch raise, lower or recline chair up to 180° without interrupting the exam. An auto-return button returns patient to the low, upright position.

● Ultra-Durable Membrane Switches

Backlit membrane switches last 10-times longer than mechanical units and run at low voltages. A programmable safety function prevents unwanted operation of the chair.

● Finishes and Upholstery

- Hand-sewn upholstery covers high-density foam that lasts for years.
- Removable cushions for easy cleaning or replacement.
- Base available in #20 Pearl.
- Choose from 14 standard vinyl colors or unlimited custom selections – call Reliance at 1-800-735-0357 or visit my color studio for a color card.

● CSA Mark of Quality

The 710 has achieved the CSA, C/US certificate of compliance to assure safety and durability.

● Tech Specs

Base: 24¼" diameter. Depth in chair position (excluding headrest): 37". Table position: 64". Seat: 19" wide by 18" deep. Distance between armrests: 19". Two base heights. High: 24⁵/₈" lowered, 35¼" raised. Low: 22¹/₈" lowered, 30⁵/₈" raised. Shipping weight: 325 lbs. Electrical: 115 volts, 60/50 Hz, 6.3 amperes. Optional electric: 230 volts, 50 Hz, 4.0 amperes.